



- Low Reagent Volume ■
- Open System ■
- Auto Cuvette Detection ■
- Auto Start Function ■
- Sample Incubation Beeper ■

MISPA CLOG
NEXT GENERATION COAGULOMETER

Authorised Distributor :

SCIENTIFIC
TRADERS

T 0771 2535351, 2225373, 2227753, 4032645



TECHNICAL SPECIFICATION

Principle

Turbodensitometric with automated start function.
Opto-mechanical principle with automatic zero adjustment and magnetic stir bar for homogenizing of the test suspension and increased sensitivity.

Measuring Unit

Single channel measuring with magnetic mixer motor.
Measuring block heated, controlled at 37.4°C +/- 0.3°C.
Single reagent position with stirring function.
4 cuvette positions.

Tests (installed)

PT
aPTT
Fibrinogen
Thrombin Time
Single Factors

Test Volume

Micro Cuvette Min. 150µl - max. 300µl.
7 pre-programmed test positions.
Automated calculation of seconds, INR, %, g/l, mg/dl, Ratio.
Calibration-/ reference curves editable up to 9 meas. points.
Automated cuvette detection.
Double or single determinations.

Display

LCD Single line with 8 characters.
Alphanumeric membrane keypad.

Power supply

AC 100V-240V, 47-63Hz, 260mA
Output: 12VDC 800mA

Weight

0.7 kg

Dimensions


190 x 130 x 60 mm (WxDxH)

Options

External Printer

MEASURING METHOD

The analyser operates according to the opto-mechanical measuring principle. This measuring principle is especially suited for lipemic and/or icteric colored samples as well as reagents with kaolin. A light beam passes through the cuvette containing the test plasma onto a photo detector. Any change in the intensity of the transmitted light, that is light increase or decrease, is converted into an electric signal. Hence, even the most unstable clot can be detected. The period from adding the start reagent until clot formation is measured. It then can be converted into the appropriate units

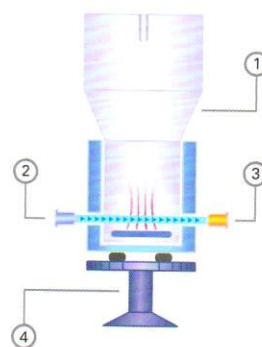

1. Test Cuvette
2. Lamp
3. Detector
4. Stirrer Motor

AGAPPE DIAGNOSTICS LTD.

Corporate Off: / Factory: 'Agappe Hills', Pattimattom (PO), Dist. Ernakulam, Kerala - 683 562, India.
TEL: + 91 484 2689992, 2687401 FAX: + 91 484 2689666 productcorp@agappe.in | www.agappe.com

The company reserves the right to change any design and technical features of the product at anytime, if needed.

AGAPPE 
Your Best Partner In Diagnostics