

ABX Micros ES 60

Hematology analyzer



PHYSICAL SPECIFICATIONS

Dimensions & Weight:

Height	Width	Depth	Weight OT	Weight CT
16.9 in	14.2 in	14.2 in	30.9 lbs	35.7 lbs
43 cm	36 cm	36 cm	14 kg	17 kg

Throughput*:

OT: 60 samples/hour

CT: 50 samples/hour

* Throughput performance is reachable in optimized conditions with an admitted margin for variation of 10%.

Sound Pressure Level:

< 53 dB

Operating Temperature & Humidity:

16 to 30°C (61 to 86°F) room temperature.

Maximum relative humidity 85% for temperature up to 30°C (86°F)

Specimen Volume:

CBC + 3 part DIFF 10 µL

Power Requirements:

Power supply from 100 V to 240 V (± 10%)
50 Hz to 60 Hz

Power consumption 100 VA (-30%, +10%)

Reagents:

- ABX Minidil LMG (10L or 20L)

- ABX Minilyse LMG (1L)

- ABX Miniclean (1L)

- ABX Minoclar (0.5L)

- ABX Cleaner (0.5L)

- Minilysebio (0,4 and 1L)

METHODS & TECHNOLOGIES

RBC/PLT Detection Principles

Method	Impedance
Ruby aperture diameter	50 µm
Depression of count	200 mb
Duration of count	2 (or 3) x 6 seconds
Dilution ratio	1/15 000

HGB Measurement

Method	Photometry
Wavelength	550 nm
Dilution ratio	1/260

HCT Measurement

Method	Numeric integration
--------	---------------------

WBC Measurement

Method	Impedance
Ruby aperture diameter	80 µm
Depression of count	200 mb
Duration of count	2 (or 3) x 6 seconds
Dilution ratio	1/260

3 Part Differentiation

Method	Impedance + specific lyse action
--------	----------------------------------

MCV, MCH, MCHC, RDW-CV, RDW-SD*, PCT*, PDW*

Calculation

SOFTWARE

8.4" LCD screen

256,000 colors, 640 x 480 pixel resolution

Inverter card for LCD screen

Touch screen technology

Secure digital memory card: 512 Mb

58mm thermal printer with rapid paper loading

Linux Operating System

PARAMETERS & PERFORMANCE DATA

19 Parameters:

WBC	RBC	MCH	PDW*
LYM# & LYM%	HGB	MCHC	PLT
MON# & MON%	HCT	RDW-SD*	MPV
GRA# & GRA%	MCV	RDW-CV	PCT*

Linearity:

Parameter	Linearity Limit	Visible Range
WBC	0 - 100 x 10 ⁹ /mm ³	100 - 150 x 10 ⁹ /mm ³
RBC	0 - 8 x 10 ⁶ /mm ³	8 - 18 x 10 ⁶ /mm ³
HGB	0 - 26 g/dL	26 - 30 g/dL
HCT	0 - 80 %	80 - 90%
PLT (Whole blood)	0 - 2200 x 10 ⁹ /mm ³	2200 - 6000 x 10 ⁹ /mm ³
PLT (PLT concentrate)	0 - 4000 x 10 ⁹ /mm ³	4000 - 6000 x 10 ⁹ /mm ³

Precision:

Parameters	%CV	Range
WBC	<2.5	4-10.0 x 10 ⁹ /µL
RBC	<2.0	4-7 x 10 ⁶ /µL
HGB	<1.5	12-18 g/dl
HCT	<2.0	36-54%
PLT	<5.0	200-500 x 10 ⁹ /µL

Carry-over:

WBC	<2%
RBC	<2%
HGB	<2%
HCT	<2%
PLT	<2%

* RUO parameters (Research Use Only)

CERTIFICATION

EN 61326-1:2006

EN 61326-2-6: 2006

IEC 61010-1: 2010

IEC 61010-2-081: 2001 + A1

IEC 61010-2-101: 2002

UL 61010-1: 2012

CAN/CSA - C22.2 N.61010-1-12

CAN/CSA - C22.2 N.61010-2-081-04

CAN/CSA - C22.2 N.61010-2-101-04

IVD: 98/79/EC

HORIBA
Medical

HORIBA India Private Limited

246, Okhla Industrial Estate Phase III, New Delhi - 110020, India, Tel: +91 (11) 4646 5000 / 4669 5001 Fax: +91 (11) 4646 5020 / 4669 5010
http://www.horiba.com, Toll Free No. 1800-103-4470, E-mail: pentramicros.hin@horiba.com

FRANCE +33 (0)4 67 14 15 15 - ITALY +39 / 06 51 59 22 1 - SPAIN +34 / 91- 353 30 10 - PORTUGAL +351 / 2 14 72 17 70 - UK +44 (0) 1604 542650
POLAND +48 / 22 30 05 460 - USA +1 / 949 453 0500 - BRAZIL +55 / 11 2923-5439 - THAILAND +66 / 2 861 59 95 - CHINA +86 / 21 3222 1818
INDIA +91 / 11 4646 5000 - GERMANY AXON LAB AG +49 / 7153 92260 - DISTRIBUTORS NETWORK +33 (0)4 67 14 15 16

HORIBA Medical online : http://www.horiba.com/medical



IMPRIM'VERT®

HORIBA
Medical

ABX Micros ES 60

The global analyzer for Indian pathologies



Advanced PAMS technology
Real time reagent Management
Comprehensive erythrocyte indices
HLM-7 compatibility
QC flexibility

Subject to modifications / Ref. SAP: 1300034557 / © 2017 HORIBA Medical - HORIBA ABX SAS - France - RCS Montpellier 328 031 042 (02/2017)

Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA

Explore the future

Automotive Test Systems | Process & Environmental | Medical | Semiconductor | Scientific

HORIBA

ABX Micros ES 60

The Gold Standard

A generation of hematology analyzer based on the Micros concept, universally recognized for its reliability, robustness and high quality results

More Powerful, More Intuitive, More Secured, More Connectivity...

The perfect combination of intuitive use and innovative technologies



- 60 samples per hour (OT Version)
- CBC + 3 DIFF (19 parameters)
- Large color touch screen
- Comprehensive quality control program
- New cyanide free lyse reagent


Open tube version
(Compatible with some closed tubes)



Integrated barcode reader*

- . Secured identification (patient, QC and reagent lot numbers)
- . Compatible with ISBT128 barcodes

*Available in closed tube version
(optional external barcode reader for open tube version)

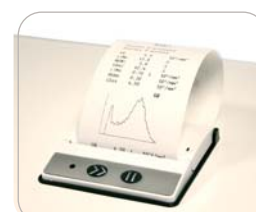

Integrated ticket printer


- . Full data print-out
- . PLT, RBC, and WBC curves

(optional external printer in closed tube and open tube versions)



USB connection

- . Archiving patients and QC results
- . Downloading QC target values

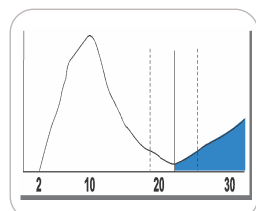


Advanced PAMS technology

- . Platelet Aggregate Management System
- . Accurate Platelet counting by distinguishing PLT and Microcytes

- User friendly interface with rapid access to all menus and virtual keyboard
- User profile management with different access levels
- Flexible data exchange for patient results and quality control
- Multi-connection capability (network and peripherals)

Easy interpretation

- Control all the parameters from the main screen
- Access to detailed results (normal ranges and curves)
- 1000 patients results capacity
- Unlimited patient report archiving (external viewer)
- Customized reporting format

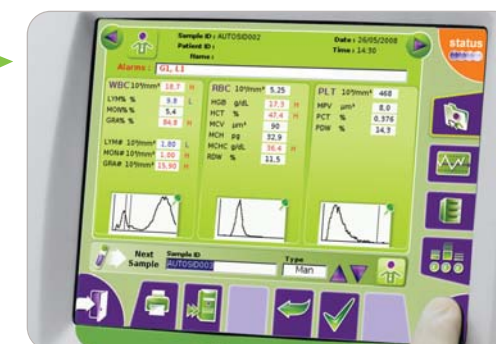

Comprehensive quality control program

- 3 simultaneous activated blood control levels
- Levey-Jennings graphics
- XB management
- Import the 3 QC levels data in one time
- Overlapping QC capabilities (2 active lots)



Automatic system monitoring

- Real-time instrument status with help messages
- Automatic reagent test count down
- Automatic daily maintenance and operational alerts
- Complete traceability in compliance with accreditation standards (users, reagents, controls, maintenance...)



ABX Micros Concept "Zero maintenance" system

- **Micro-sampling**
Only 10µl of whole blood per analysis
- **Liquid valves**
Precise volume and reliability
- **Stepper motor vacuum pump**
Noiseless running (no-compressor)

Full Connectivity Integrated Technologies

- **Software**
Java™, Linux Operating System, HTML
- **Hardware: 6 connections ports**
RS232 (2), USB (3), RJ45 (1)
- **ABX, ASTM and HL7 protocols available**

