

www.sysmex-ap.com

Sysmex Corporation

1-5-1, Wakino-hama-Kaigandori, Chuo-ku, Kobe,
Hyogo 651-0073, Japan
Tel +81 (78) 265-0500 Fax +81 (78) 265-0524
www.sysmex.co.jp

Sysmex Asia Pacific Pte Ltd

9, Tampines Grande #06-18,
Singapore 528735
Tel +65 6221-3629 Fax +65 6221-3687
www.sysmex-ap.com

Sysmex India Pvt Ltd

308 Ascot Centre, 3rd Floor, Sahar Airport Road,
Andheri (East), Mumbai-400099, India
Tel +91 (22) 2822-4040 Fax +91 (22) 2836-5068
www.sysmex.co.in

PT Sysmex Indonesia

Cyber 2 Tower, 5th Floor, Unit E,
Jl. HR. Rasuna Said Blok X5 No. 13,
Jakarta Selatan 12950, Indonesia
Tel +62 (21) 3002-6688 Fax +62 (21) 3002-6699
www.sysmex.co.id

Sysmex (Malaysia) Sdn Bhd

11A & 15, Jalan PJS 7/21, Bandar Sunway,
46150 Petaling Jaya, Selangor, Malaysia
Tel +60 (3) 5637-1788 Fax +60 (3) 5637-1688
www.sysmex.com.my

Sysmex New Zealand Ltd

382 - 386 Manukau Rd, PO Box 26-085, Epsom,
Auckland 1344, New Zealand
Tel +64 (9) 630-3554 Fax +64 (9) 630-8135
www.sysmex.co.nz

Sysmex (Thailand) Co., Ltd.

12th Floor, Tonson Tower, 900 Ploenchit Road, Lumpini,
Pathumwan, Bangkok 10330, Thailand
Tel +66 (2) 635-5777 Fax +66 (2) 257-0530
www.sysmex.co.th

Sysmex Philippines Inc

30th Floor MDC 100 Bldg, E. Rodrigues Jr. Ave. cor.
Eastwood Ave. Bagumbayan, Quezon City, Philippines
Tel +63 (2) 621-2460 Fax +63 (2) 621-2432
www.sysmex.com.ph

Sysmex Vietnam Co. Ltd.

Unit 802, Centre Point Building, 106 Nguyen Van Troi St,
Phu Nhuan Dist, Ho Chi Minh City, Vietnam
Tel +84 (8) 3997-9400 Fax +84 (8) 3997-9405
www.sysmex.com.vn



We Believe the Possibilities.

Simplifying coagulation testing

CA-101

Semi-automated blood coagulation analyzer

About Sysmex

Sysmex is a global market leader in the development and implementation of clinical diagnostic and health IT products and services for laboratories, hospitals and healthcare organizations.

We deliver total solutions in the field of clinical laboratory testing, including hematology, hemostasis and urinalysis. Supplying products and services to customers in more than 150 countries, Sysmex is the market leader in hematology and coagulation instrumentation worldwide.

By working together with our customers and business partners, Sysmex is committed to our mission of shaping the advancement of healthcare.



CA-101 Your Straightforward, Compact Coagulation Analyzer



Set cuvette



Dispense and Incubate sample



Timer countdown



Aspirate start reagent



Dispense start reagent, measurement is started automatically. Result is displayed on the LCD screen or printed out automatically to an external printer

Ideal for laboratories that require minimum manual intervention in performing routine coagulation assays. It is designed to simplify coagulation testing; making it more standardized and reproducible

- Coagulation testing is made easy with built-in timer, temperature control, and auto-start and timing of measuring process
- Easy to use with pre-programmed methods
- Objective detection of clot formation
- Automatic light intensity adjustment according to the turbidity of the plasma, making it possible to measure icteric or lipemic plasma reliably
- Sophisticated standard curve tool to maintain accuracy and precision
- Increase throughput with short incubation time
- A technology that requires only half the sample and reagent volume for testing as compared to manual method
- A compact design, combined with the use of ChipCARD system for your convenience



ChipCARD

The analyzer allows you to load method-specific data and parameters for a certain method via an integrated ChipCARD reader


Turbodensitometric Measuring Principle

A light beam passes through the cuvette containing the test plasma onto a photo detector. Any change in the intensity of the transmitted light, i.e., light increase or decrease, is converted into an electric signal

The stirrer mix the reagent and plasma in the cuvette. And at the same time creates a small whirl through the mixer movement to ensure even the smallest fibrin clot is formed in front of the photo detector

Once the start reagent is added, the lamp intensity automatically adjusts up or down according to the turbidity of the sample

This stirring action and optical measurement constitute the basic features of the turbodensitometric measuring principle

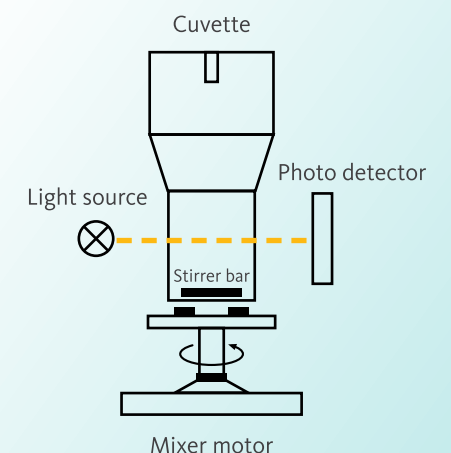

Together with high quality reagents, it provides a reliable coagulation solution, not only for low volume testing, but also as an backup for fully automated system that requires a different optical method

Reagent system

- | | |
|-------------|--|
| PT | • Thromborel S, Dade Innovin |
| APTT | • Dade Actin, Dade Actin FS, Dade Actin FSL, Pathromtin SL |
| TT | • Thromboclotin, Test Thrombin Reagent |
| Fbg | • Multifibren U |